start[®]
by iHealth

START BY IHEALTH THF
Non contact infrared thermometer
Model: MDI901

Instruction for use
Instructions d'utilisation
Benutzerhandbuch
Instrucciones de uso
Istruzioni per l'uso
Instruções de utilização
Gebruiksaanwijzing

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Thank you for choosing our Start by iHealth THf non-contact infrared thermometer. This device is used to measure an object's temperature based on the relationship between temperature and measurable infrared radiation. Simply aim the unit's probe toward the surface to be measured to obtain a quick and accurate temperature. To ensure proper use, please be sure to read this user manual carefully, paying close attention to the safety precautions.

SAFETY PRECAUTIONS

- Use of this thermometer is not intended as a substitute for consultation with your physician. It is dangerous for user to perform a self-evaluation and self-treatment based on the measuring result. Be sure to follow doctor's instruction.
- Keep the thermometer out of reach of children. For accidental swallow of battery or other component, please consult the doctor at once.
- Don't throw the battery into fire.

NOTICE

- The device is precision instrument, don't drop, tramp or impose any vibration or impact on the thermometer.
- Do not touch the lens of the probe with your fingers, and do not disassemble the device by yourself.
- Before measuring forehead temperature, make sure the hair re-moved, sweat dried.
- After you do some exercise, eating and bathing, you should stay still indoor about 30 minutes before measurement.
- To make the measurement data reliable and stable, when ambient temperature varies a lot, the thermometer should be placed indoors for about 30 minutes before using.
- When we measure somebody continuously, the temperature should be measured every minute, if you need to measure yourself continuously for a short time, there are some slight errors when you read the temperature ,which is a normal phenomenon. At this time, we should choose the average. We recommend that you measure yourself no more than 3 times in a specific unit of time, because the temperature of the human will conduct to the thermometer, it may affect the accuracy of measurement.
- There is no absolute standard about the temperature of the human, so please try to collect the recording of individual temperature in the usual, as a reference for having a fever or not.

- Do not measure the sites of scarred tissue or tissue compromised by skin disorders, because sensing body temperature from sites of scarred tissue or tissue compromised by skin disorders.
- Do not measure the site of forehead temperature if that patient has trauma on forehead.
- Do not measure if that patient is treated with certain drug therapies.
- Do not immerse the device into water or any other liquid, and not directly sunlight exposure.
- Do not use a mobile or cordless phone near the thermometer when measuring.
- Body temperature may increase in the drug within the effort time limit. Please don't measure.
- In order to ensure the accuracy of measurement data, please don't take measurement of body temperature in strong electromagnetic interference environment (such as microwave, high frequency equipment operation environment).
- Do not disassemble, repair, or modify the unit.
- Do not touch any metal part in the battery compartment during measuring.
- The thermometer must be stored according to the technical specifications.
- The materials (ABS) of expect contact with patient had passed the ISO 10993-5 and ISO 10993-10 standard test, no toxicity, allergy and irritation reaction. They are compliance with the MDD requirements. Based on the current science and technology, other potential allergic reactions are unknown.
- The patient is an intended operator. The patient can measure, read data and replace battery under normal circumstances and maintain the device and its accessories according to the user manual.
- The device is not intended for PATIENT transport outside a healthcare facility.

RECOMMENDATIONS

- Don't use this thermometer for other purposes.
- It is forbidden to leave the product exposed to any chemical solvent, direct sunshine or high temperature.
- Don't expose the thermometer under direct sunlight long time so as not to damage the battery.
- Do not measure while talking on the phone.
- Please report to MANUFACTURER if any unexpected operation or events occur.

INTENDED USE

This thermometer is intended to measure forehead temperature at home or hospital, including anyone, such as infants, children and adults.
For the safety reason, children or the baby's temperature must be measured by parent or adults.

MEASUREMENT MODES

- Forehead temperature.

Measure the skin surface of human forehead's temperature accurately, take the place of traditional mercury thermometer and electrical thermometer.

Normal forehead temperature range for different ages

Ages (years)	Normal temp. (°C)	Normal temp. (°F)
0-2	36.4-38.0	97.5-100.4
3-10	36.1-37.8	97.0-100.0
11-65	35.9-37.6	96.6-99.7
>65	35.8-35.5	96.4-99.5

NOTE: The normal temperature and difference between the different body parts is individual. To define yours, measure your temperature for at least 2 weeks at the same forehead position and time.

NOTE: When consulting your physician, communicate that which position the temperature is measured by infrared thermometer, note the individual's normal temperature range by infrared thermometer as additional reference.

NOTE: Because the forehead temperature is affected obviously by the external environment (eg: environment, air convection and skin tone, etc), we advice that you take the forehead temperature only as reference.

- Object temperature.

Measure the surface temperature of object or liquid, such as bathwater, milk temperature etc.

- Room temperature.

BACKLIGHT COLOR INDICATION



When the correct result is displayed, the backlight from LCD will light different color (6 seconds) to indicate the temperature range as below:

Temp. Range	Color
$T < 37.5^{\circ}\text{C}$ ($T < 99.5^{\circ}\text{F}$)	Green
$37.5 \leq T \leq 38.0^{\circ}\text{C}$	Orange
$38.0^{\circ}\text{C} \leq T$	Red

MEMORY

The thermometer stores the last 99 temperature measurements. The last temperature will be automatically displayed when it is turned on again. The display will show "M" symbol.

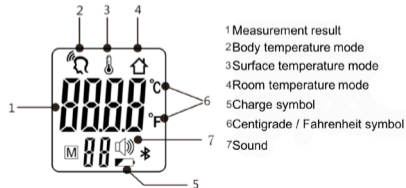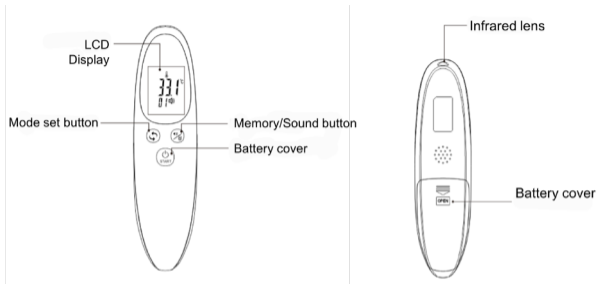DESCRIPTION

- 1 Measurement result
- 2 Body temperature mode
- 3 Surface temperature mode
- 4 Room temperature mode
- 5 Charge symbol
- 6 Centigrade / Fahrenheit symbol
- 7 Sound

MODE 	Press for mode selection / Previous memory check
MEM / SOUND 	Check memory data / Next memory check / Sound switch
START 	Measuring key / Power ON / Confirming key

OPERATION INSTRUCTION

PREPARATION

1) Check sensor

If have pollution and spray, please clean it. (The cleaning method see the chapter 8 Care and Cleaning for detailed.) If the sensor's lens is damaged, please stop using.

2) Check thermometer

When you press "Start" button, the system will have self testing of software and hardware. If there are problems, LCD will display "Err" symbol. If found, please refer to the chapter Trouble-shooting.

3) In order to make the accurate measuring result, put the thermometer in the measurement environment for 30 minutes.

4) Unexpected fluctuations in ambient temperature may decrease the accuracy of the measurement results. When the thermometer at the same measurement position display at different ambient temperature, or measure temperature in front of the air conditioner, it will not be able to obtain accurate results.

5) If you measure forehead temperature, clean forehead and arrange hair, make sure the forehead is naked and clean, in order to ensure the accuracy of measurement.

BATTERY LOADING

The thermometer is supplied with two AAA batteries. Insert new batteries when the battery symbol appears on the display. Open the battery compartment. Remove the batteries and replace with new batteries, making sure the poles are in the right direction. Slide battery door back until it snaps in place.

Notice:

- Please observe the related national laws of disposing the abandoned battery and don't litter to the garbage can.
- Please take out the batteries if the device is not used for long periods of time.

- Please don't put the batteries in the fire.

To protect the environment, dispose of empty battery at your retail store or at appropriate collection sites according to national or local regulations.

MODE CHOICE

Press "Start" button to turn on. Press "Mode" button to toggle between the three modes.

BODY TEMPERATURE MEASUREMENT

- Choose the Body temperature mode
- Aim the forehead, press "Start" button, the temperature is displayed immediately.

Take the measurement as shown:

To ensure accurate readings please refer to the guidelines below.

*When taking body temperatures, the optical focus should be aimed at the center of the forehead and 1-3cm away from the forehead. Hold the thermometer in a vertical position. Ensure the forehead is free from hair.

*If sweating has occurred, avoid taking forehead temperatures. In these circumstances, it is recommended that temperatures are taken from the armpit.

*Take at least three readings and take an average.

*If the thermometer is taken from an extremely cold place to an extremely hot place or vice versa, it is recommended to wait 30 minutes to ensure accurate measurements.

*Do not use under intensive sunlight.

*Wait 15 minutes after being outdoors in cold or hot weather before taking temperature readings.

*If your baby is crying, wait for the crying to cease before taking a reading.

* If the environmental temperature varies significantly, place the thermometer in the room for 15-20 minutes before taking temperature.

* Do not expose the thermometer to direct sunlight, or place in wet or damp conditions for prolonged periods.

SURFACE TEMPERATURE MEASUREMENT

- Choose the surface temperature mode
- Aim the object surface, press "Start" button, the temperature is displayed immediately.

ROOM TEMPERATURE MEASUREMENT

- Choose the room temperature mode
- Aim the room, press "Start" button, the temperature is displayed immediately.

MUTE / TURN ON THE SOUND

When thermometer is off, press "MEM/Sound" button to choose "OFF" or "ON" to mute or turn on

the beep sound.

DISPLAYING STORED RESULTS

Press "Start" button to turn on the thermometer, press the "MEM/Sound" button to activate the recording function. Press "MEM/Sound" button again to review previously recorded temperature readings.

CARE AND CLEANING

The probe tip and lens are the most delicate part of the thermometer. It has to be clean and intact to ensure accurate readings.

If the thermometer is ever accidentally used, please clean the probe and lens as follows:

- Very gently wipe the surface with a cotton swab or soft cloth moistened with alcohol. Only use the thermometer after the alcohol has completely dried out.
- If the lens is damaged, contact the distributor.

Clean the unit body:

- Use a soft, dry cloth to clean the thermometer display and unit body.
- If very dirty, use a soft cloth with alcohol to cleaning.

NOTES:

- Don't use abrasive cleaners.
- Don't use other non-recommended methods to perform disinfect.
- Non-waterproof, don't use the abrasive cleaner to clean the product, don't drop the thermometer in the water or the other liquid.

MAINTENANCE

- We do not authorize any institution or individual to maintain and repair of the product. If you suspect there is something wrong on the products, please contact the manufacturer or distributor.
- The user must not attempt any repairs to the device or any of its accessories. Please contact the retailer for repair.
- Opening of the equipment by unauthorized agencies is not allowed and will terminate any claim to warranty.

CALIBRATION

The thermometer is initially calibrated at the time of manufacture. If this thermometer is used

according to the use instruction, periodic re-adjustment is not required. If any time your question the accuracy of measurement, please contact distributor or manufacturer.

STORAGE

- Don't put the thermometer under the sunshine, high temperature and moist environment or some-place which maybe get in touch with fire or is vulnerable to vibration.
- Take out the battery if don't use the device in a long time.

TROUBLESHOOTING

Message	Possible cause	Solution
"Hi"	Temp. higher than 43°C (109.4°F).	Please leave the status and wait for 30 minutes to measure. Re-measure according to the manual.
"Lo"	Temp. lower than 32.0°C (89.6°F).	
"ErH"	The ambient temperature is higher of range of measurement (40°C/104°F).	Keep the thermometer in the room where temperature is (10°C-40°C/50°F -104°F) for 30 minutes.
"ErL"	The ambient temperature is lower of range of measurement (10°C/50°F).	
"Err"	The hardware is damaged or thermometer system fails or affected by electric magnetic field.	Excluding the possibility of temperature allowance first , then send the device to your dealer for repair.
	Low battery.	Replace the batteries.

SPECIFICATIONS

- Model name: MDI901
- Size: 161x42x30mm
- Weight: 79g (without batteries)
- Power supply: 2x 1,5VDC AAA
- Battery life: ~1000 measurements
- Measurement accuracy: $\pm 0.2^{\circ}\text{C}/0.4^{\circ}\text{F}$ between 35.5~42.0°C (95.0~107.6°F);
-5- $\pm 0.3^{\circ}\text{C}/0.5^{\circ}\text{F}$ between 32.0~35.4°C (86.9~95.7°F) and 42.1~43°C (107.8~109.4°F)

- Repeated testing deviation: $\leq 0.3^{\circ}\text{C}$ (0.5°F)
- Testing distance: 1-3cm
- Testing time: 1-2s
- Measuring range:
 - Body temp. mode: $32^{\circ}\text{C}\sim 43^{\circ}\text{C}$ ($89.6^{\circ}\text{F}\sim 109.4^{\circ}\text{F}$)
 - Surface temp. mode: $0^{\circ}\text{C}\sim 100^{\circ}\text{C}$ ($32^{\circ}\text{F}\sim 212^{\circ}\text{F}$)
- Operating conditions:
 - $10\sim 40^{\circ}\text{C}$ ($50\sim 104^{\circ}\text{F}$),
 - humidity $\leq \text{Rh}95\%$, 70-106KPa
- Storage conditions:
 - $-25\sim 55^{\circ}\text{C}$ ($-13\sim 131.0^{\circ}\text{F}$),
 - humidity $\leq \text{Rh}95\%$, 70-106KPa

PACK CONTENTS

- 1x Start by iHealth THf
- 1x User manual
- 2x AAA batteries

STANDARDS

This unit fulfils the provisions of the EC directive 93/42/EEC (Medical Device Directive) and the European Standards.

EN ISO 80601-2-56:2012: Particular requirements for basic safety and essential performance of clinical thermometers for body temperature measurement.

ELECTROMAGNETIC COMPATIBILITY

This product has been tested and approved according to EN60601-1-2-2007 on electromagnetic compatibility. But we cannot guarantee that the product will not be interfered by electromagnet. Do not use in a strong electromagnetic environment.

WARRANTY INFORMATION

iHealth Labs Europe ("iHealth") warrants the Start by iHealth THf electronic thermometer (the "Product"), and only the Product, against defects in materials and workmanship under normal use for

a period of two years from the date of purchase by the original purchaser ("Warranty Period"). Under this Limited Warranty, if a defect arises and a valid claim is received by iHealth within the Warranty Period regarding the Product, at its option and to the extent permitted by law, your retailer will either (1) repair the Product using new or refurbished replacement parts or (2) exchange the Product with a new or refurbished Product. In the event of a defect, to the extent permitted by law, these are the sole and exclusive remedies.

Start by iHealth is a trademark of iHealthLabs Europe.

iHealth is a trademark of iHealthLabs Europe.

This warranty does not apply: (a) to consumable parts, such as the battery that diminish over time, unless failure has occurred due to a defect in materials or workmanship; (b) to cosmetic damage, including but not limited to scratches, dents ; (c) to damage caused by accident, abuse, misuse, contact with liquid; (d) to damage caused by operating the Start by iHealth product outside the user manual, the technical specifications or other Start by iHealth product published guidelines; (e) to damage caused by service performed by anyone who is not a representative of iHealth or one of its representatives

To contact the technical support:

- European support:

support@ihealthlabs.eu

If you have questions or need assistance outside the operational days and times, please contact your health care provider

Shenzhen Medek Bio-Medical Co., Ltd.

No.101 Wenhao Industrial Park, No.13 Yuanxin Road, Tongle Community, Baolong Street, Longgang District, 518100 Shenzhen, P.R.China

Tel: +86-755-61996096

Fax: +86-755-61996096

Website: www.medekchina.com

Shanghai International Holding Corp. GmbH (Europe)

Eiffestrasse 80, 20537 Hamburg, Germany

Tel: +49-40-2513175

Fax: +49-40-255726

E-mail: shholding@hotmail.com

EXPLANATION OF SYMBOLS

 Authorized representative in the European Community

 Serial number

 Consult instructions for use

 Manufacturer

 Environmental Protection–Electrical products waste should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for

recycling advice.

 Keep in a dry place

 Caution! Consult accompanying documents

 Batch code

 Storage temperature Limit

 Type BF applied parts

 Complies with the requirements of the European Directive 93/42/EC

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ShenZhen Medek Bio-Medical Co., Ltd.